

# Wound drainage systems

NSC has a range of high quality wound drainage systems that can be used for all types of procedures and surgeon's preferences

- Privac-Low: low vacuum wound drainage systems
- Privac-Tvac: low vacuum thumbhole wound drainage systems
- Privac-High: high vacuum wound drainage systems
- Privac-Mini: small vacuum wound drainage systems
- Privac-Reinfusion autologous blood transfusion systems
- Surgiblake silicone channel drains

ISO 9001  
BUREAU VERITAS  
Certification



National Surgical Corporation  
is certified to  
AS/NZS ISO 9001:2008 (NSW & WA)



BE SUPPORTED BY AN  
AUSTRALIAN COMPANY  
PROVIDING VALUE BEYOND THE PRODUCT

**NATIONAL SURGICAL CORPORATION**

*A proudly founded and owned Australian company - providing value to our customers beyond the product*

# Wound drainage systems

*National Surgical Corporation have a range high quality wound drainage systems: cost effective, low or high suction, closed systems, silicone drains*

## **VARIOUS SUCTION OPTIONS**

NSC wound drains are available in low or high suction which combines gravity drainage to active drainage

## **IDEAL FOR GENERAL OPERATIONS**

Ideal for general operations as well as specialised procedures including orthopaedic, mastectomies, gynaecology and plastic surgery

## **NON RETURN VALVE TUBING ON LOW VACUUM**

The NSC low vacuum drains feature a non-return valve on the bellows to ensure fluids move away from the wound and drain into bag

## **COMPLETE KITS AVAILABLE**

All NSC wound drains are available in complete kits including bellows/reservoir, universal connection, 10FG, 14FG, 18FG drain tube and trocar

## **CLOSED SYSTEM**

Complies with Australian Infection Control Standards and OH&S guidelines

## **STERILE REPLACEMENT BAGS/RESERVOIRS**

Sterile, single use replacement bags and reservoirs are available to use post-operatively on the ward

## **PATIENT COMFORT & MOBILITY**

Hanger allows attachment of system close to patient, on clothing or mobility aid

## **SILICONE CHANNEL DRAIN**

The NSC Surgiblake silicone channel drain is flexible and soft which ensures less tissue trauma on removal of drain, making it more comfortable for the patient

## **COMPATIBILITY**

The Surgiblake silicone drain can be used with any low vacuum drain including Privac-Low & Privac-Tvac

## **WOUND DRAINAGE OPTIONS**

The NSC range includes Privac-Low, Privac-Tvac, Privac-High, Privac-Mini and Privac-Reinfusion wound drainage systems and Surgiblake silicone channel drains

## **PRIVAC-LOW**

- Low suction: combining gravity drainage to a maximum of 118mm hg/ 16 kPa active drainage
- Non-return valve on 275ml bellows to ensure fluids move away from the wound and drain into bag
- Use with PVC or silicone drains

## **PRIVAC-TVAC**

- Low suction: combining gravity drainage to a maximum of 90mm hg/ 12 kPa active drainage
- Non-return valve on 275ml bellows to ensure fluids move away from the wound and drain into bag
- Use with PVC or silicone drains
- Bellows with thumbhole for quick and easy vacuuming with one hand



## PRIVAC-HIGH

- High suction: maximum suction to 700mmHg/ 93 kPa
- For use in any surgery where a high vacuum drain is preferred. Up to 600ml capacity
- Use with PVC drains
- Different reservoir sizes: 200ml, 400ml, 600ml



## PRIVAC-MINI

- Low suction: maximum suction to 150mm Hg/ 20 kPa
- Ideal for minor operations requiring less fluid drainage including head, neck, hand or plastic surgery
- 50ml capacity bellows, easy to squeeze/vacuum




## PRIVAC-REINFUSION

- Low suction: maximum suction to 120mm hg/ 16 kPa
- Reduces the need for homologous blood transfusions by enabling the patient to have their own filtered blood reinfused
- Control the flow rate during reinfusion with 40µm blood filter, integrated dip chamber and rollerclamp



## SURIGBLAKE SILICONE CHANNEL DRAINS

- Manufactured with a single piece extrusion which reduces risk of breakage compared with two piece silicone drain designs
- Flexible and soft silicone ensures less tissue trauma on removal of drain, making it more comfortable for the patient
- Four flow channels for superior tissue contact, no holes reduces tissue ingrowth and less risk of blockage
- Available as a round or flat drain, with or without trocar and connector to suit most wound drain systems
- Can be used with any low vacuum system including Privac-Low or Privac-Tvac wound drainage systems

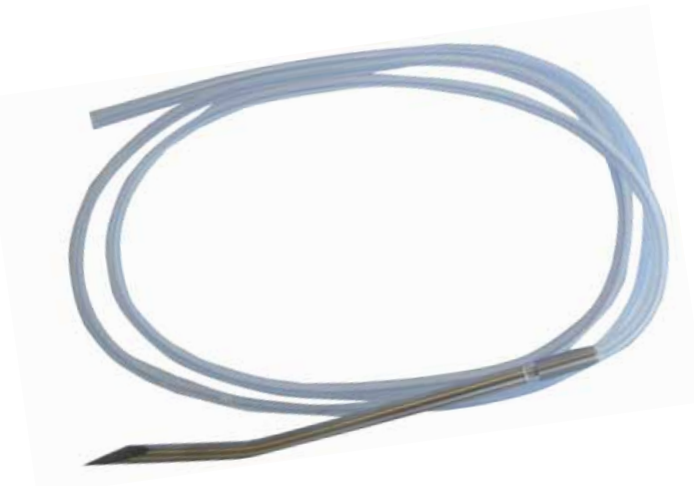

# Let NSC support you



**BE SUPPORTED BY AN  
AUSTRALIAN COMPANY**

PROVIDING VALUE BEYOND THE PRODUCT

## What makes NSC different to our competitors to support you?

- A proudly Australian owned and founded company
- Combining 25 years of expertise in the healthcare industry
- Creating value for our customers beyond the product
- Working with our customers to offer specialised solutions
- Offering a large product range focused on cost effectiveness, patient care and productivity
- Low overheads allows us to provide cost savings to our customers compared to the multinationals
- Coordinated, whole company approach to every customer
- Establishing trust for long term customer relationships

## Who are we?

National Surgical Corporation (NSC) is an import manufacturer in the healthcare industry, distributing surgical instruments, consumables, equipment and shelving all major hospitals and healthcare facilities throughout Australia. Our extensive product range can be used by all specialist disciplines and departments within a hospital, day surgery or specialist centre.

## Let NSC support you and your department

Don't be afraid to change, NSC will support and assist you throughout every step. NSC has an extensive training and inservicing material available and can tailor a program to meet your individual needs so that all personnel are aware of how the product operates. Switching to NSC can save you significantly and we can show you how!

## Our range includes



Wound drains



Electrosurgical  
range



OT and general  
consumables



Surgical  
instruments



Laparoscopic  
instruments



Medichill  
ice packs



Protection  
products



Cleaning &  
sterilisation



Shelving &  
storage



High density  
storage

## **CUSTOMER SERVICE**

**Phone: 1800 138 138 (free call)**

**Fax: 1800 138 140 (free fax)**

**[info@nationalsurgical.com.au](mailto:info@nationalsurgical.com.au)**

**[www.nationalsurgical.com.au](http://www.nationalsurgical.com.au)**

**Contact your National Surgical Corporation Account Manager for more information**



NATIONAL SURGICAL CORPORATION

Customer Service 1800 138 138 • [www.nationalsurgical.com.au](http://www.nationalsurgical.com.au) • [info@nationalsurgical.com.au](mailto:info@nationalsurgical.com.au)